



Medical Technologies

Easy on-PC

Portable & PC spirometer


Connecting patients to precise pulmonary care

ndd's PC-driven spirometer employs the interactive **EasyOne Connect** software to deliver rapid, accurate and reliable results. With compelling features such as live curve visualization, animated incentives and integration with leading EMR/EHR systems, the Easy on-PC is the clear choice for PC-connected spirometry.

ndd's patented **TrueFlow** technology eliminates the necessity for calibration, saving time while producing precise and dependable results.




Simple. Certain. Proven.

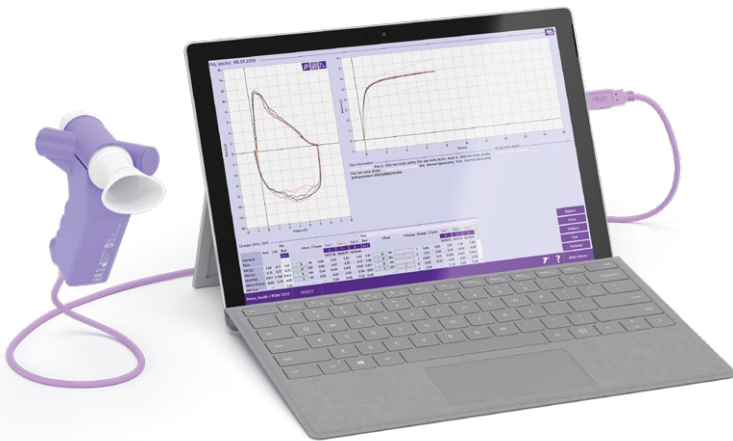

Incentive screens

Available with **EasyOne Connect** software

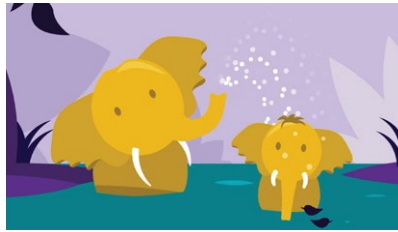

Specifications

Test Types: FVC, FVL, SVC, MVV, Challenge and Pre/Post testing

Languages: Chinese, Croatian, Danish, Dutch, English, Finnish, French, German, Italian, Japanese, Norwegian, Portuguese, Russian, Spanish, Swedish, Turkish, Vietnamese

EMR connectivity: HL7, XML, GDT, API with EasyOne Connect

Standards: Meets 2005 & 2019 ATS/ERS standards and NIOSH/OSHA

Hard disk capacity: Installation/system 1 GB. Data up to 4 GB

RAM: 2 GB

EasyOne Connect software compatibility: Windows 7, Windows 8 and 8.1 (32 and 64 bit), Windows 10 (32 and 64 Bit), Windows 11

Article numbers

2700-3: Easy on-PC spirometer

2050-1: Spirette box of 50

2050-5: Spirette box of 200

2050-10: Spirette box of 500

Features

- Real-time animated incentives
- Proven integration with top EMR/EHR systems
- Long-term trending of results
- Automated quality feedback and interpretation
- Intuitive software interface
- Integrated provocation test with customizable protocols
- Wide range of selectable predicted values
- Calibration-free technology

TrueFlow
makes the difference

nnd's unique ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity.

nnd **TrueFlow** is a resistance-free solution that does not require calibration during its lifetime.

EasyOne Connect

EasyOne Connect software offers one database and one platform for all **EasyOne** point-of-care solutions. Data management has never been easier.



nnd Medical Technologies
300 Brickstone Square, Suite 604 Andover, MA, 01810 USA
Phone: +1-978-470-0923

nnd Medizintechnik AG
Technoparkstrasse 1, CH-8005 Zurich, Switzerland
Phone: +41 44 512 65 00



SCAN ME

