




# EasyOne Pro/Pro LAB

## Complete point of care pulmonary function testing

Portable DLCO, lung volumes & spirometry.  
A highly efficient and accurate way to test lung function.

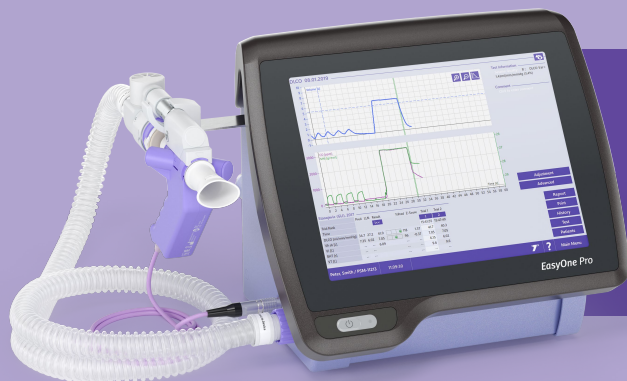

### EasyOne Pro

- Spirometry
- DLCO
- Provocation testing

#### Lung volumes

- Single breath helium dilution

# Simple. Certain. Proven.

## Features

- Automated user guidance & test quality feedback based on 2019 & 2005 ATS/ERS standards
- Z-score, LLN and % predicted for fast interpretation of results
- Real-time curves and pediatric incentives
- Export of pdf files and raw data
- Flexible HL7 and XML interface for easy EMR integration
- Hygienic solution with Spirette and Barriette consumables eliminates the risk of cross-contamination
- Compact device with smooth surfaces for easy and thorough cleaning

### EasyOne Pro LAB

- Spirometry
- DLCO
- Provocation testing

#### Lung volumes

- Multi-breath nitrogen washout



## Specifications

**EasyOne Pro:** FVC, FVL, SVC, MVV, Pre/Post, Provocation, DLCO

**EasyOne Pro LAB:** FVC, FVL SVC, MVV, Pre/Post, Provocation, DLCO, MBW, LCI

**Languages:** English, Chinese, Croatian, Danish, Dutch, Finnish, French, German, Italian, Japanese, Norwegian, Portuguese, Russian, Spanish, Swedish, Turkish, Vietnamese

**EMR connectivity:** HL7, XML, GDT, API with EasyOne Connect

**Dimensions:** 10.6 x 13.2 x 10.6" / 17.65 lb / 27 x 33.5 x 27 cm (H x W x D), 8 kg

## Part numbers

**3000-1:** EasyOne Pro

**3100-1:** EasyOne Pro LAB

**2050-1:** Spirette box of 50

**2050-5:** Spirette box of 200

**2050-10:** Spirette box of 500

**3050-1:** DLCO Barriettes box of 50

**3050-2:** DLCO Barriettes box of 100

**3150-1:** FRC Barriettes box of 40 (for use with EasyOne Pro LAB only)

**3150-2:** FRC Barriettes box of 80 (for use with EasyOne Pro LAB only)

**TrueFlow**  
makes the difference

The original ultrasonic flow measurement is highly accurate in all flow ranges, independent of gas composition, pressure, temperature and humidity.

ndd **TrueFlow** is a hygienic and resistance-free solution that does not require calibration during its lifetime.

## EasyOne Connect

**EasyOne Connect** software offers one database and one platform for all EasyOne point-of-care-solutions. Data management has never been easier.

**TrueCheck**  
automated precision

**TrueCheck** takes care of the essential quality control for gas analysis testing. EasyOne Pro is the only device proven to be accurate for a lifetime for DLCO measurements.



SCAN ME

